
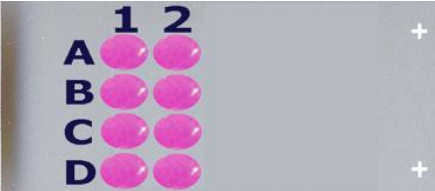


Cat No: NDRG1081 – NDRG1 IHC control Array

Lot#	Cores	Size	Cut	Format	QA/QC	 	
NDRG1081	8	1.5mm	4um	4X2	H&E, IHC anti-NDRG1		

Recommended applications: NDRG1 IHC control with strong, moderate, low/negative expressers.

Description: The array contains eight 1.5mm cores duplicated from four cases of appropriate tissue samples. They were fixed in 10% neutral buffered formalin for 24 hours and processed using identical SOPs. Sections were picked onto Superfrost Plus or APES coated Superfrost slides. They can be stored for use at 4C for up to six months from the date of shipment. **There may be 5 to 10% of tissue core loss.** As validated by IHC, A1, 2 is a strong (>3) expresser, B1, 2 is a moderate to strong (2~3) expresser, C1, 2 is weak to moderate (1~2) expresser while D1, 2 is strong and non-expresser of the NDRG1 molecule.

Array position	Sex	Age	Anatomic site	Pathology	Nature	Grade	Stage	NDRG1 Staining
A01	F	58	Uterus, cervix	Squamous cell carcinoma	Malignant	III	T1N0M0	3.3~3.5
A02	F	58	Uterus, cervix	Squamous cell carcinoma	Malignant	III	T1N0M0	3.3~3.5
B01	F	31	Uterus, cervix	Squamous cell carcinoma	Malignant	II~III	T2N0M0	2.0~2.5
B02	F	31	Uterus, cervix	Squamous cell carcinoma	Malignant	II~III	T2N0M0	2.0~2.5
C01	F	41	Uterus, cervix	Squamous cell carcinoma	Malignant	I~II	T1N1M0	1.0~2.0
C02	F	41	Uterus, cervix	Squamous cell carcinoma	Malignant	I~II	T1N1M0	1.0~2.0
D01	M	38	Skin	Normal	Normal			Derma 0, Epidermis 3.0
D02	M	38	Skin	Normal	Normal			Derma 0, Epidermis 3.0

Notes: Bake at 60C for ~60 minutes before use. If antigen retrieving is needed, it is important to avoid **direct-boiling and high pH or high strength** antigen retrieving buffer. For availability of complimentary IHC data, please contact us at info@pantomics.com.

Certified by: Langxing Pan, M.D., Ph.D.

Interpretation of Pantomics IHC Analysis

"0" is negative. "0.5" is borderline staining with no significance. "1" is weak staining. "1.5" is weak staining with foci of moderate staining. "2" is moderate staining. "2.5" is moderate staining with foci of strong staining. "3" is homogeneous strong staining. "3.5" is very strong and homogeneous staining with no significant background. "4" is over staining usually with background staining.