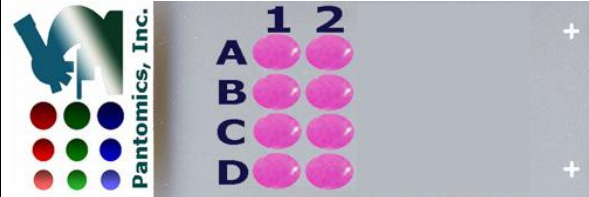


**Cat No: TFF3081 – Trefoil factor 3 (TFF3) IHC control Array**

Lot#	Cores	Size	Cut	Format	QA/QC		
TFF3081	8	1.5mm	4um	2X4	H&E, IHC anti-Trefoil factor 3		

**Recommended applications:** Trefoil factor 3 (TFF3) IHC control with strong, moderate, low/negative expressers.

**Description:** The array contains eight 1.5mm cores duplicated from four cases of appropriate tissue samples. They were fixed in 10% neutral buffered formalin for 24 hours and processed using identical SOPs. Sections were picked onto Superfrost Plus or APES coated Superfrost slides. They can be stored for use at 4C for up to six months from the date of shipment. **There may be 5 to 10% of tissue core loss.** As validated by IHC, A1, 2 is a strong (>3) expresser, B1, 2 is a moderate to strong (2~3) expresser, C1, 2 is weak to moderate (1~2) expresser while D1, 2 is weak or non-expresser of the Trefoil factor 3 (TFF3) molecule.

Array position	Sex	Age	Anatomic site	Pathology	Nature	Grade	Stage	Trefoil factor 3 Staining
A01	F	41	Breast	Invasive ductal carcinoma	Malignant	I~II	T2N0M0	3.2~3.5
A02	F	41	Breast	Invasive ductal carcinoma	Malignant	I~II	T2N0M0	3.2~3.5
B01	M	66	Intestine, colon	Adenocarcinoma	Malignant	II	T3N0M0	2.5~3.0
B02	M	66	Intestine, colon	Adenocarcinoma	Malignant	II	T3N0M0	2.5~3.0
C01	M	66	Intestine, colon	Adenocarcinoma	Malignant	II	T3N0M0	1.0~2.0
C02	M	66	Intestine, colon	Adenocarcinoma	Malignant	II	T3N0M0	1.0~2.0
D01	M	83	Prostate	Hyperplasia	Benign			0~0.5
D02	M	83	Prostate	Hyperplasia	Benign			0~0.5

**Notes:** Bake at 60C for ~60 minutes before use. If antigen retrieving is needed, it is important to avoid **direct-boiling and high pH or high strength** antigen retrieving buffer. For availability of complimentary IHC data, please contact us at [info@pantomics.com](mailto:info@pantomics.com).

**Certified by:** Langxing Pan, M.D., Ph.D.

### Interpretation of Pantomics IHC Analysis

"0" is negative. "0.5" is borderline staining with no significance. "1" is weak staining. "1.5" is weak staining with foci of moderate staining. "2" is moderate staining. "2.5" is moderate staining with foci of strong staining. "3" is homogeneous strong staining. "3.5" is very strong and homogeneous staining with no significant background. "4" is over staining usually with background staining.