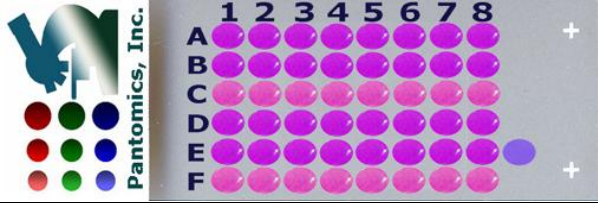


Cat No: COC481 – Colon cancer/adjacent uninvolved tissue array

Lot#	Cores	Size	Cut	Format	QA/QC	
COC48101	48	2mm	4um	6X8	H&E, IHC anti-Cytokeratin	

Recommended applications: For Research use only. Designed for colon normal vs cancer RNA or protein tissue profiling using IHC or ISH. Duplicated tumor tissues from each case are arrayed to minimize heterogeneity.

Description: 16 cases of colon cancer (purple), each in duplicates, with corresponding uninvolved tissue (pink) from the same patient as control.

All the tissues were from surgical resection. They were fixed in 10% neutral buffered formalin for 24 hours and processed using identical SOPs. Sections were picked onto Superfrost Plus or APES coated Superfrost slides. They all have a guaranteed six months' shelf-life at 4C from the date of shipment. There may be >5% tissue core loss, but the core retaining rate should be >90%.

Position	Sex	Age	Anatomic Site	Pathology	Grade	Stage
A01	F	60	Intestine, colon	Adenocarcinoma	I~II	T2N1M0
A02	F	85	Intestine, colon	Adenocarcinoma	I~II	T2N0M0
A03	F	43	Intestine, colon	Adenocarcinoma	I~II	T2N0M0
A04	M	51	Intestine, colon	Adenocarcinoma	II	T2N0M0
A05	M	65	Intestine, colon	Adenocarcinoma	II	T3N0M0
A06	F	48	Intestine, colon	Adenocarcinoma	II~III	T3N0M0
A07	F	74	Intestine, colon	Adenocarcinoma	II	T2N1M0
A08	F	43	Intestine, colon	Adenocarcinoma	II~III	T4N1M0
B01	F	60	Intestine, colon	Adenocarcinoma	I~II	T2N1M0
B02	F	85	Intestine, colon	Adenocarcinoma	I~II	T2N0M0
B03	F	43	Intestine, colon	Adenocarcinoma	I~II	T2N0M0
B04	M	51	Intestine, colon	Adenocarcinoma	II	T2N0M0
B05	M	65	Intestine, colon	Adenocarcinoma	II	T3N0M0
B06	F	48	Intestine, colon	Adenocarcinoma	II~III	T3N0M0
B07	F	74	Intestine, colon	Adenocarcinoma	II	T2N1M0
B08	F	43	Intestine, colon	Adenocarcinoma	II~III	T4N1M0
C01	F	60	Intestine, colon	Uninvolved colon tissue of A1,B1	null	null
C02	F	85	Intestine, colon	Uninvolved colon tissue of A2,B2	null	null
C03	F	43	Intestine, colon	Uninvolved colon tissue of A3,B3	null	null
C04	M	51	Intestine, colon	Uninvolved colon tissue of A4,B4	null	null
C05	M	65	Intestine, colon	Uninvolved colon tissue of A5,B5	null	null
C06	F	48	Intestine, colon	Uninvolved colon tissue of A6,B6	null	null
C07	F	74	Intestine, colon	Uninvolved colon tissue of A7,B7	null	null
C08	F	43	Intestine, colon	Uninvolved colon tissue of A8,B8	null	null
D01	M	36	Intestine, colon	Adenocarcinoma	III	T3N0M0
D02	M	47	Intestine, colon	Adenocarcinoma	III	T3N1M0
D03	F	53	Intestine, colon	Adenocarcinoma	III	T3N1M0
D04	F	53	Intestine, colon	Adenocarcinoma	III	T3N2M0
D05	M	72	Intestine, colon	Adenocarcinoma	III	T2N0M0
D06	M	36	Intestine, colon	Adenocarcinoma	III~IV	T2N0M0
D07	F	60	Intestine, colon	Adenocarcinoma	II~III	T3N0M0
D08	M	35	Intestine, colon	Adenocarcinoma	III~IV	T4N2M1
E01	M	36	Intestine, colon	Adenocarcinoma	III	T3N0M0
E02	M	47	Intestine, colon	Adenocarcinoma	III	T3N1M0
E03	F	53	Intestine, colon	Adenocarcinoma	III	T3N1M0
E04	F	53	Intestine, colon	Adenocarcinoma	III	T3N2M0
E05	M	72	Intestine, colon	Adenocarcinoma	III	T2N0M0
E06	M	36	Intestine, colon	Adenocarcinoma	III~IV	T2N0M0



Advancing Biomedical Science Through Tissue Arrays

E07	F	60	Intestine, colon	Adenocarcinoma	II~III	T3N0M0
E08	M	35	Intestine, colon	Adenocarcinoma	III~IV	T4N2M1
F01	M	36	Intestine, colon	Uninvolved colon tissue of D1,E1	null	null
F02	M	47	Intestine, colon	Uninvolved colon tissue of D2,E2	null	null
F03	F	53	Intestine, colon	Uninvolved colon tissue of D3,E3	null	null
F04	F	53	Intestine, colon	Uninvolved colon tissue of D4,E4	null	null
F05	M	72	Intestine, colon	Uninvolved colon tissue of D5,E5	null	null
F06	M	36	Intestine, colon	Uninvolved colon tissue of D6,E6	null	null
F07	F	60	Intestine, colon	Uninvolved colon tissue of D7,E7	null	null
F08	M	35	Intestine, colon	Uninvolved colon tissue of D8,E8	null	null

Notes: Bake at 60C for 30 minutes before use. If antigen retrieving is needed, it is always a good idea to start with a protocol with weak to mild strength. Please check our TNM guide below as a reference. For IHC related information on products that provide them, please contact us at info@pantomics.com for the most up-to-date information

Certified by: Langxing Pan, M.D.

TNM Classification: Colon and Rectum Carcinoma

T- Primary tumor

TX - Primary tumor cannot be assessed

TO - No evidence of primary tumor

Tis - Carcinoma in situ; intraepithelial or invasion of lamina propria

T1 - Tumor invades submucosa

T2 - Tumor invades muscularis propria

T3 - Tumor invades through muscularis propria into subserosa or into non-peritonealized pericolic or perirectal tissues.

T4 - Tumor directly invades other organs or structures and/or perforate visceral peritoneum

N - Regional lymph nodes

NX - Regional lymph nodes cannot be assessed

NO - No regional lymph node metastasis

N1 - Metastasis in 1 to 3 regional lymph nodes

N2 - metastasis in 4 or more regional lymph nodes

M - Distant metastasis

MX - Distant metastasis cannot be assessed

MO - No distant metastasis

M1 - Distant metastasis