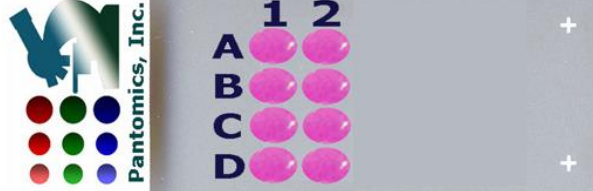


Cat No: TYROSIN081 – Tyrosinase IHC control Array

Lot#	Cores	Size	Cut	Format	QA/QC		
TYROSIN081	8	1.5mm	4um	4X2	H&E, IHC anti-Tyrosinase		

Recommended applications: Tyrosinase IHC control with strong, moderate, low/negative expressers.

Description: The array contains eight 1.5mm cores duplicated from four cases of appropriate tissue samples. They were fixed in 10% neutral buffered formalin for 24 hours and processed using identical SOPs. Sections were picked onto Superfrost Plus or APES coated Superfrost slides. They can be stored for use at 4C for up to six months from the date of shipment. **There may be 5 to 10% of tissue core loss.** As validated by IHC, A1, 2 is a strong (>3) expresser, B1, 2 is a moderate to strong (2~3) expresser, C1, 2 is weak to moderate (1~2) expresser while D1, 2 is weak or non-expresser of the Tyrosinase molecule.

Array position	Sex	Age	Anatomic site	Pathology	Nature	Grade	Stage	Tyrosinase Staining
A01	F	54	Intestine, rectum	Melanoma	Malignant		T2N0M0	3.3~3.5
A02	F	54	Intestine, rectum	Melanoma	Malignant		T2N0M0	3.3~3.5
B01	M	45	Lymph node, mesenteric	Metastatic Melanoma	Metastasis			2.5~3.0
B02	M	45	Lymph node, mesenteric	Metastatic Melanoma	Metastasis			2.5~3.0
C01	F	35	Skin	Nevus	Benign			1.5~2.0
C02	F	35	Skin	Nevus	Benign			1.5~2.0
D01	M	58	Skin	Squamous cell carcinoma	Malignant	I	T1N0M0	0~0.5
D02	M	58	Skin	Squamous cell carcinoma	Malignant	I	T1N0M0	0~0.5

Notes: Bake at 60C for ~60 minutes before use. If antigen retrieving is needed, it is important to avoid **direct-boiling and high pH or high strength** antigen retrieving buffer. For availability of complimentary IHC data, please contact us at info@pantomics.com.

Certified by: Langxing Pan, M.D., Ph.D.

Interpretation of Pantomics IHC Analysis

"0" is negative. "0.5" is borderline staining with no significance. "1" is weak staining. "1.5" is weak staining with foci of moderate staining. "2" is moderate staining. "2.5" is moderate staining with foci of strong staining. "3" is homogeneous strong staining. "3.5" is very strong and homogeneous staining with no significant background. "4" is over staining usually with background staining.